




Olympic Brainz Monitor
Exporting Data in CSV or EDF+ Format

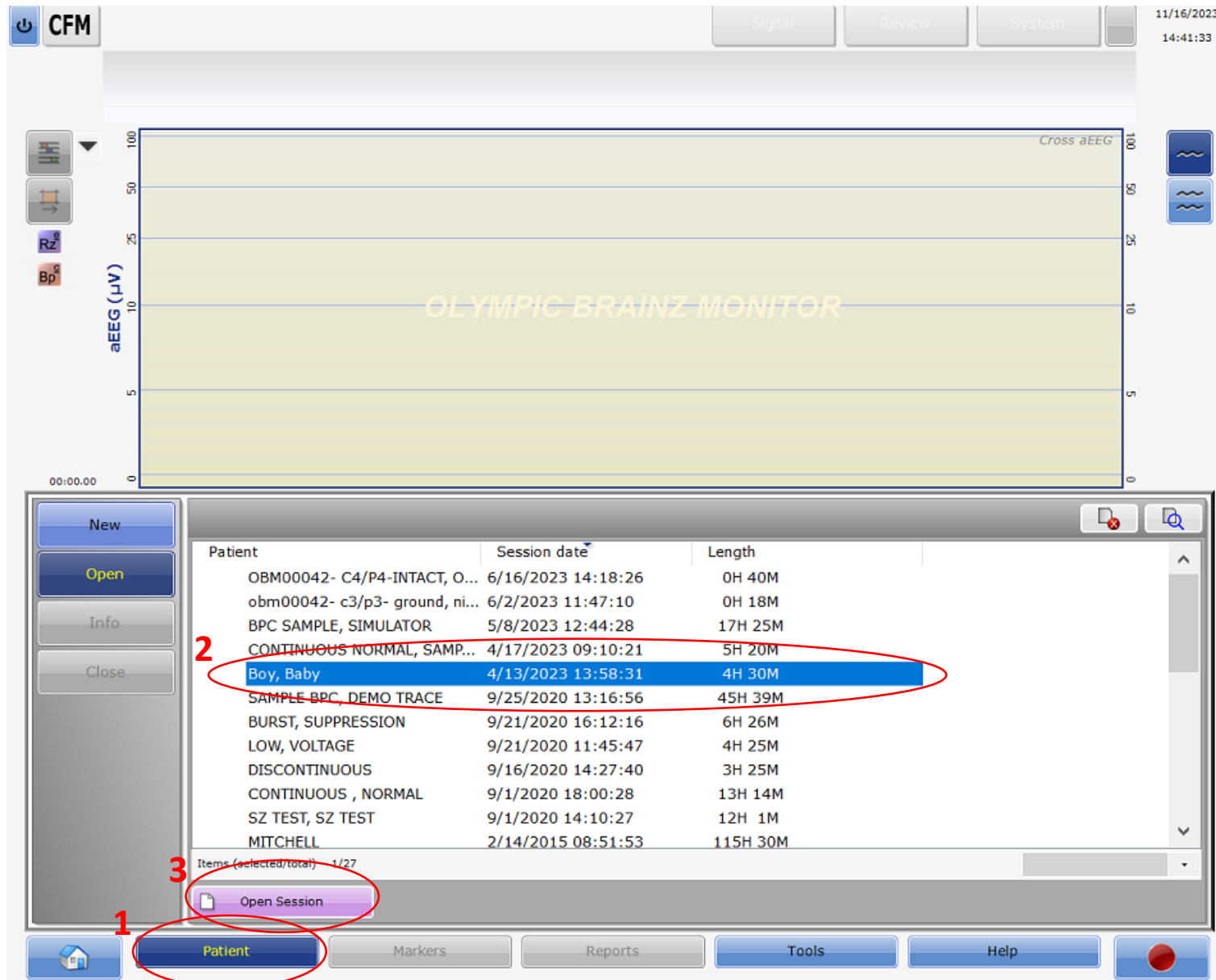


1. From the home screen, select the PATIENT function button
2. Select the patient
3. Touch on the OPEN SESSION button

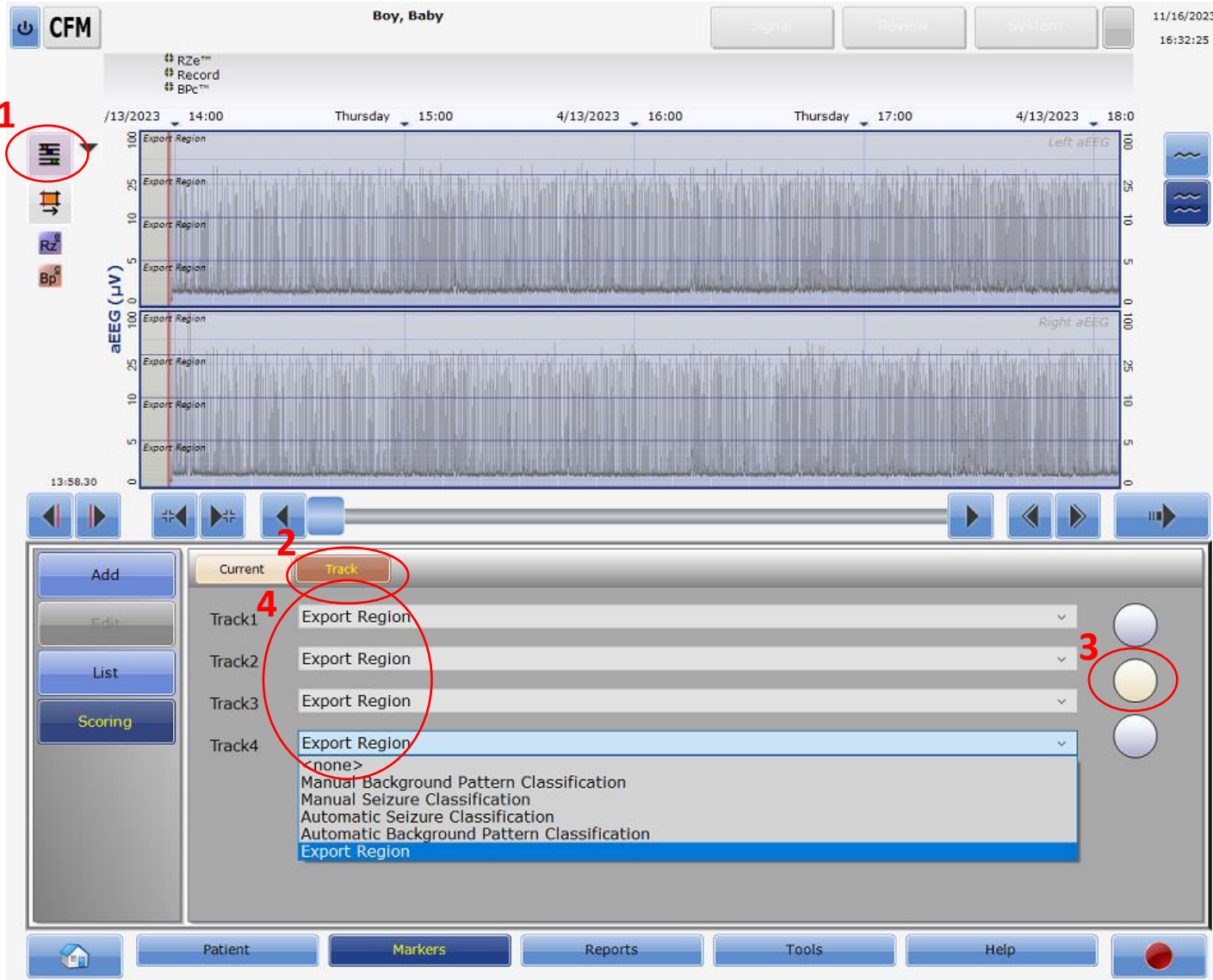


1. Select the SCORING icon
2. In the Scoring window, select TRACK
3. There are three track setting options, select the second or third option by touching on the circle
4. Set tracks 1-4 to EXPORT REGION

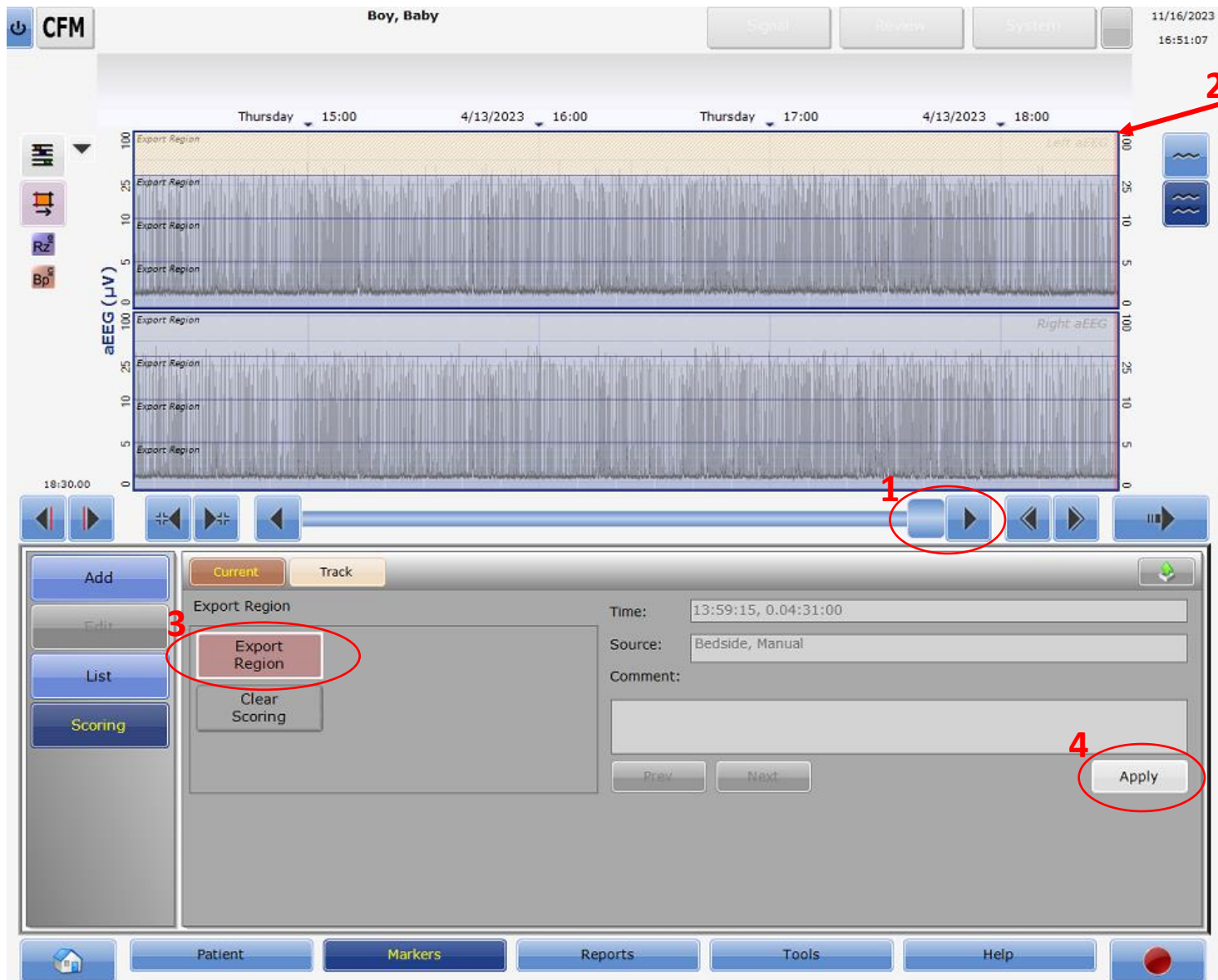


1. Drag the SCROLL BAR to the right
2. Put the cursor on the end of the recording- Will turn the track orange
3. Select the EXPORT REGION button
4. Then touch on APPLY- This will turn all EXPORT REGION tracks red

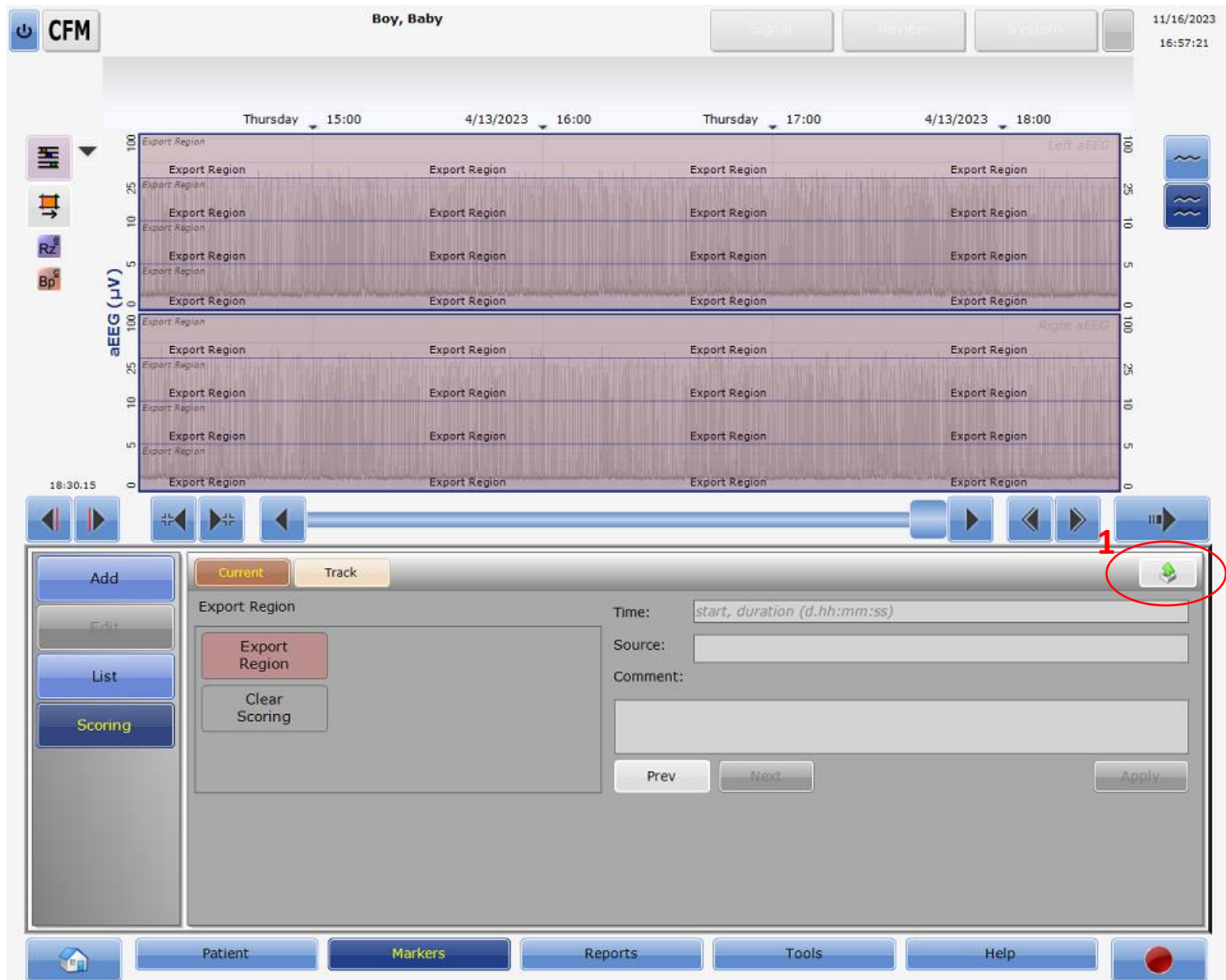

1. Select the EXPORT icon

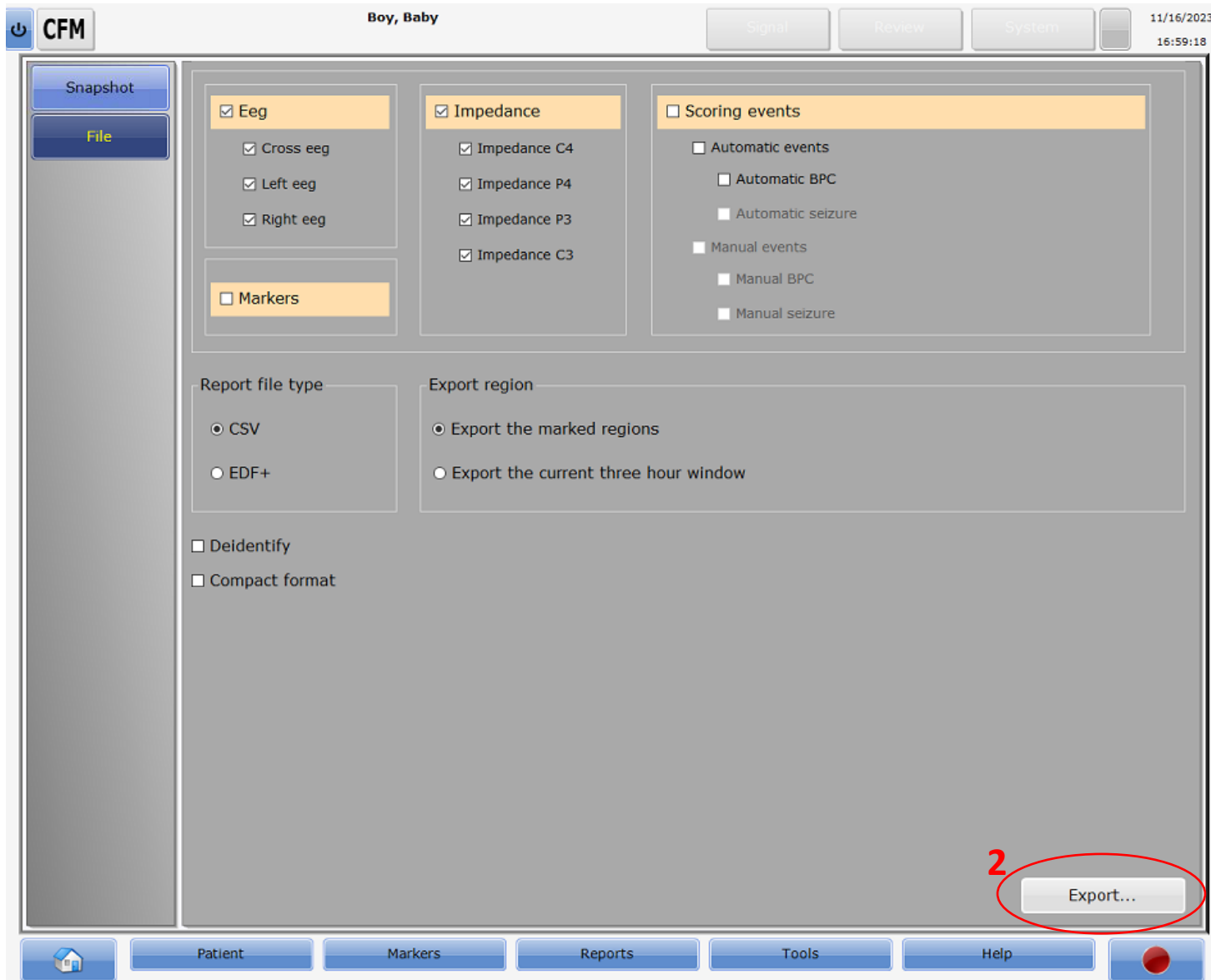

1. Check all desired data to be exported
2. Select the EXPORT icon

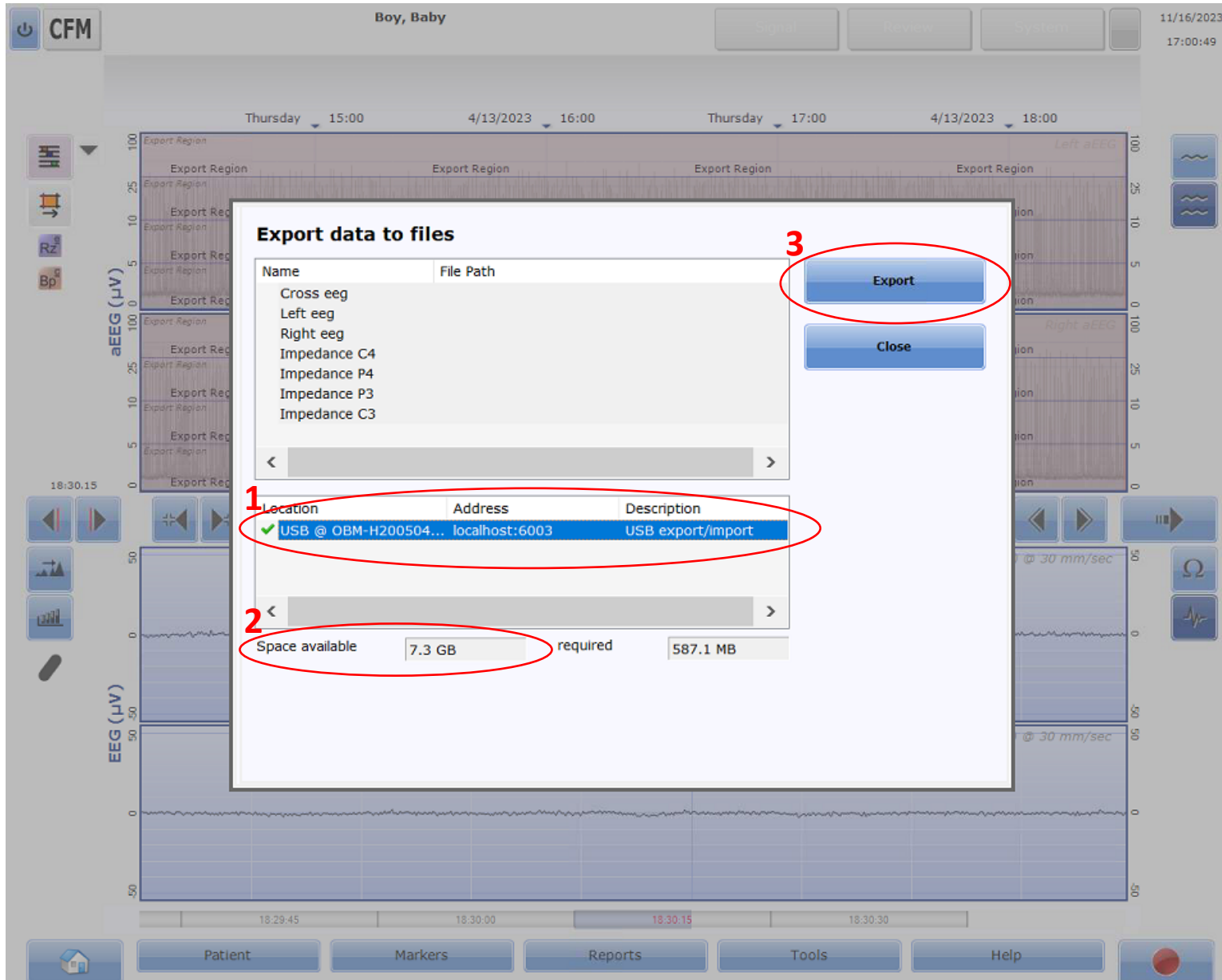

1. Select the USB drive for transfer
2. Verify available space on the drive
3. Touch on EXPORT

natus®